


ENT CHAIRS

Supplied with head-rest

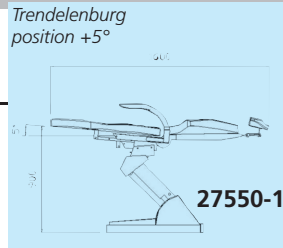

Backrest synchronized to the leg-rest movement can be adjusted at 90° and inclined until 190°

The seat is adjustable (42 cm vertical excursion)

Width 650 mm

Revolving and removable arm-rests in light grey expanded polyurethane

Foot switch for up and down movements

CHAIRS AVAILABLE IN 36 COLOURS

*See colour at page 125 Available in 60 days

Optional footrest on request

27551

Seat base in heavy metal painted in epoxydic powders, with support of pantograph arm - 85° left and right rotation - lever locking device

1.160 mm
490 mm
650 mm

27546

- 27545 OTO/PV ENT CHAIR with head support - colour on request*
- 27546 OTO/PV ENT CHAIR with head support - green Melbourne
- 27547 OTO/PV ENT CHAIR with head support - sky blue Parigi


Patient chair with Peltex-lined seat and backrest. Electric up and down movements by silent engine controlled through foot switch. Back-rest can be adjusted up to 180° by a gas spring. Chair rotates 360° and is provided with a rotation locking device placed at the bottom. Power: 230 V - 50 Hz (Available 110 V - 60 Hz). Made in Italy.

- 27550 OTOPEX ENT CHAIR with head support - colour on request*
- 27551 OTOPEX ENT CHAIR with head support - blue Chicago
- 27552 OTOPEX ENT CHAIR with head support - green Toronto

ENT chair designed also for ambulatory operations. Moved by very silent electrical engines operated by joystick foot-control. Also provided with button to re-set to previously selected working position. Upholstery made of simulated leather, high density undeformable foam padding (50 kg/cm³). Seat and leg-rest are joined in a single piece. Power: 230 V - 50 Hz (Available 110 V - 60 Hz). Made in Italy.

ENT CABINETS

Support pole for monitor

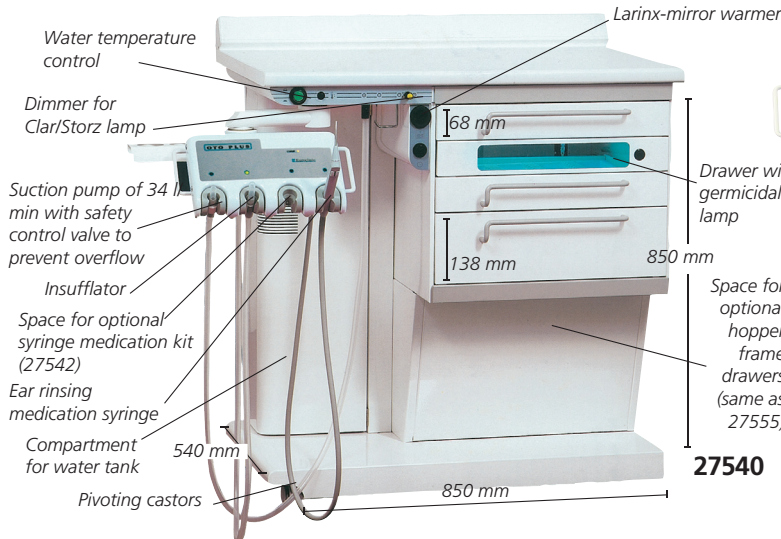


- 27540 OTOPLUS - ENT UNIT

The E.N.T. unit has rubber pivoting castors, a working surface made of scratch-proof rolled sections, 2 small and one large drawers provided of removable internal partitions, one drawer with inside germicidal lamps to keep instrument sterilization. Otoplus is provided with a silent suction pump of 34 l/min, with power supply complete with dim adjustment device for Clar/Storz lamp, larynx-mirror warmer, twin-cylinder dry compressor for tubal insufflation, 2-litres distilled water tank with acoustic signal for reserve level. The articulated arm with instrument holding tray is complete with terminal suction cannula and terminal insufflation cannula, ear rinsing syringe with automatic water temperature control, (a microprocessor keeps temperature at a constant level of 37 °C). Power 230 V - 50 Hz (Available 110 V - 60 Hz). Made in Italy.

- 27539 OTOPLUS COMPLETE

Same as code 27540 with included vasister drawers (code 27541) and endoscope support (code 27543).

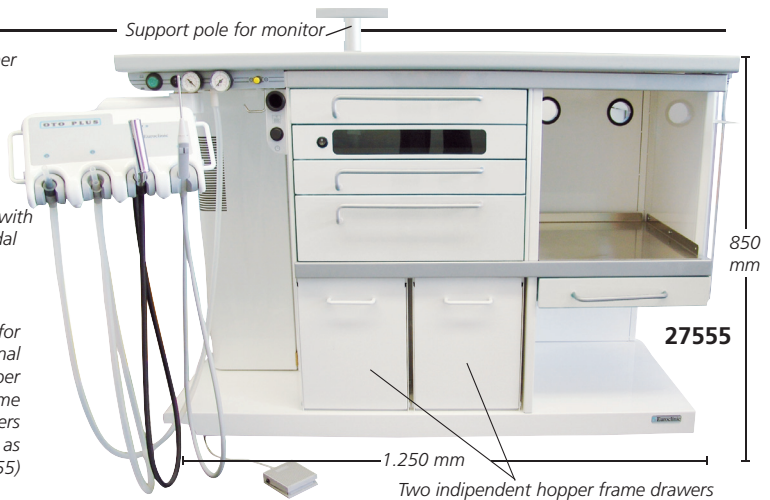

- 27554 OTOPLUS DC ENT UNIT

Otoplus unit with additional column, to be used for diagnostic equipment (light source, camera or printer) and a support pole for monitor.

- 27555 OTOPLUS DC ENT UNIT - complete

Same as code 27554 with included hopper frame drawers (code 27541) and endoscope support (code 27543)

- 27541 HOPPER FRAME DRAWERS - 2 independent drawers
- 27543 ENDOSCOPE SUPPORT - 2 transparent plastic tubes 1 tube 50 cm long and 1 tube 25 cm long both closed at the bottom.
- 27542 SYRINGE MEDICATION KIT composed of 1 syringe, 2 atomizers (1 for liquid and 1 for powder) and 3 ampoules in dark glass. To be ordered with ENT unit.
- 27544 LED LIGHT SOURCE

Makes cold white light (5500 °K, 200 lumen) for endoscopy use. It can be connected to endoscopes without need of F.O. cable. Handle connect directly to rigid endoscopes (Storz connection).

LCD display: 12.8x6.4 cm

Size: 30x25x12 cm

27544