

ERKA VARIO DESK WALL AND MOBILE - MADE IN GERMANY

- 32671 VARIO - desk - adult cuff - 285.20.481
 - 32672 VARIO - wall - adult cuff - 288.20.481+295.00.003
 - 32673 VARIO - trolley Ø 50 cm - adult cuff+288.20.481+995.00.003
- Superior quality product, made in Germany, for use in consulting rooms or hospital. The Vario model consists of a high precision manometer in a high quality ABS housing with an integrated fastening spigot for individual variable swivelling, a big cuff basket in white powder coated sheet steel, a bulb holder and spiral tubing with bulb and air release valve. This high quality chromium-plated valve makes sensitive adjustments easy. The 150 mm scale is equipped with a mirrored exterior ring and non-reflective manometer glass.
- 32898 PEDIATRIC "SUPERB" CUFF 2 TUBES - 20.5-28 cm (arm circumference)
 - 32899 ADULT "SUPERB" CUFF 2 TUBES - 27-35 cm (arm circumference) - spare
 - 32900 LARGE "SUPERB" CUFF 2 TUBES - 34-43 cm (arm circumference)

TECHNICAL SPECIFICATIONS

Height with basket:	260 mm (32672-3)
Width with basket:	180 mm (32672-3)
Depth with basket:	120 mm (32672-3)
Graduation:	0-300 mmHg ±3 mmHg
Weight with basket:	1,350 g
Height of trolley:	830-1,360 mm (adjustable)
Tubing system Ø:	18-25 mm
Swivel range:	up to 60 deg°

32671



32672



32673

ERKA.

RIESTER BIG BEN® - MADE IN GERMANY

- 32742 RIESTER BIG BEN® - rail (clamp code 32766 and rail not included)
- 32743 RIESTER BIG BEN® - wall
- Rotation 130°, 65° right, 65° left
- 32746 RIESTER BIG BEN® - trolley Ø 54 cm

Aneroid sphygmomanometer, made of resistant plastic, with large dial (147 mm) and white pointer on black dial plate for easy-to-read measurements from 0 mmHg to 300 mmHg. Basket on back for cuff storage after use. Velcro-closure black cuff 2 tubes, and long coil tubing. Latex free.

Quality product: made in Germany.

- 32766 RAIL CLAMP
- Universal clamp (for code 32742) for any rail.

Made in Germany



32743


CE 0124

Riester

Plastic basket on back for cuff storage



BHS VALIDATION

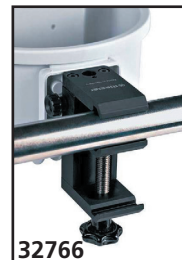

32766



CE 0124

New trolley Ø 54 cm

OMRON PROFESSIONAL B.P. MONITOR

- 32936 OMRON DIGITAL AUTOMATIC B.P. MONITOR HEM-907

Professional B.P. monitor designed for clinical use in professional settings.

- gold standard accuracy
- validated according to international protocols:
- 1- Association for the Advancement of Medical Instrumentation (AAMI)
- 2- European Society of Hypertension
- International Protocol (ESH - IP)
- intellisense technology
- irregular pulse indicator
- 3 modes of operation: single mode, average mode, auscultation (manual mode).




Supplied with battery pack, adult cuff (22-32 cm), AC adapter.

Made in Japan.

- 32937 ADULT LARGE CUFF 32-42 cm
- 32939 TROLLEY FOR OMRON HEM-907

TECHNICAL SPECIFICATIONS

Method:	Oscillometric
Measurement range:	pressure: 0-299 mmHg ± 3 mmHg or 2% pulse rate: 30-199 beats/min ± 5% of reading
Inflation:	automatic inflation with pumping
Deflation:	automatic deflation by electromagnetic control valve
Air release:	automatic rapid air release by electromagnetic control valve
Pressure detection:	electrostatic capacity semi-conductor pressure sensor
Power supply:	AC adapter (230 VAC, 50 Hz, 20 VA) Battery pack (4.8 VDC, 6W)
Electric shock protection method:	Class II B type
Dimension:	139x203x131 mm (5.1/2"x8"x5.1/6")
Weight:	approx. 910 g
Complies with:	EN 1060, EN60601-1-1 - 2007

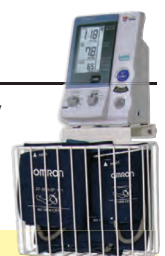

32936

+ 32939

CE 0197